


pFusionRed-N vector


The vector sequence has been compiled using the information from sequence databases, published literature, and other sources, together with partial sequences obtained by Evrogen. This vector has not been completely sequenced.

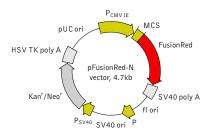

For vector sequence, please visit our Web site at http://www.evrogen.com/products/vectors.shtml

Product	Cat.#	Size
pFusionRed-N vector	FP412	20 μg
Vector type	mammalian expression vector	
Reporter	FusionRed	
Reporter codon usage	mammalian	
Promoter for FusionRed	P _{CMV IE}	
Host cells	mammalian	
Selection	prokaryotic - kanamycin eukaryotic - neomycin (G418)	
Replication	prokaryotic - pUC ori eukaryotic - SV40 ori	
Use	FusionRed expression in mammalian cells; generation of fusions to the FusionRed N-terminus	

Multiple cloning site (MCS)

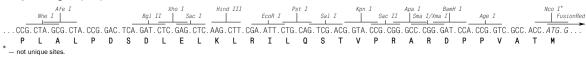

Location of features

P_{CMV IE}: 1-589 Enhancer region: 59-465 TATA box: 554-560 Transcription start point: 583 MCS: 591-671 FusionRed

Kozak consensus translation initiation site: 672-682 Start codon (ATG): 679-681; Stop codon: 1375-1377 SV40 early mRNA polyadenylation signal

Polyadenylation signals: 1530-1535 & 1559-1564 mRNA 3' ends: 1568 & 1580

f1 single-strand DNA origin: 1627-2082

Bacterial promoter for expression of Kan^r gene

-35 region: 2144-2149; -10 region: 2167-2172 Transcription start point: 2179 SV40 origin of replication: 2423-2558

SV40 early promoter

Enhancer (72-bp tandem repeats): 2256-2327 & 2328 2399

21-bp repeats: 2403-2423, 2424-2444 & 2446-2466 Early promoter element: 2479-2485

Major transcription start points: 2475, 2513, 2519 & 2524

Kanamycin/neomycin resistance gene

Neomycin phosphotransferase coding sequences: Start codon (ATG): 2607-2609; Stop codon: 3399-3401 G->A mutation to remove Pst I site: 2789

C->A (Arg to Ser) mutation to remove BssH II site: 3135 Herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signal

Polyadenylation signals: 3637-3642 & 3650-3655 pUC plasmid replication origin: 3986-4629

References

Gorman, C. (1985). "High efficiency gene transfer into mammalian cells." In: *DNA cloning: A Practical Approach, Vol. II.* Ed. by Glover. (IRL Press, Oxford, U.K.) Pp. 143–190.

Haas, J. et al. (1996) "Codon usage limitation in the expression of HIV-1 envelope glycoprotein." Curr Biol, 6 (3): 315-324 / pmid: 8805248

Kozak, M. (1987) "An analysis of 5'-noncoding sequences from 699 vertebrate messenger RNAs." Nucleic Acids Res, 15 (20): 8125–8148 / pmid: 3313277

Vector description

pFusionRed-N is a mammalian expression vector encoding red fluorescent protein FusionRed. The vector allows generation of fusions to the FusionRed N-terminus and expression of FusionRed fusions or FusionRed alone in eukaryotic (mammalian) cells.

FusionRed codon usage is optimized for high expression in mammalian cells (humanized) [Haas et al. 1996]. To increase mRNA translation efficiency, Kozak consensus translation initiation site is generated upstream of the FusionRed coding sequence [Kozak 1987]. Multiple cloning site (MCS) is located between $P_{\text{CMV}\,\text{IE}}$ and FusionRed coding sequence.

The vector backbone contains immediate early promoter of cytomegalovirus ($P_{\text{CMV IE}}$) for protein expression, SV40 origin for replication in mammalian cells expressing SV40 T-antigen, pUC origin of replication for propagation in *E. coli*, and f1 origin for single-stranded DNA production. SV40 polyadenylation signals (SV40 poly A) direct proper processing of the 3'-end of the reporter mRNA.

SV40 early promoter (P_{SV40}) provides neomycin resistance gene (Neo^r) expression to select stably transfected eukaryotic cells using G418. Bacterial promoter (P) provides kanamycin resistance gene expression (Kan^r) in *E. coli*. Kan^r/Neo^r gene is linked with herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signals

Generation of FusionRed fusion proteins

A localization signal or a gene of interest can be cloned into MCS of the vector. It will be expressed as a fusion to the FusionRed N-terminus when inserted in the same reading frame as FusionRed and no in-frame stop codons are present. The inserted sequence should contain an initiating ATG codon. FusionRed-tagged fusions retain fluorescent properties of the native protein allowing fusion localization *in vivo*. Unmodified vector will express FusionRed when transfected into eukaryotic (mammalian) cells.

Note: The plasmid DNA was isolated from dam⁺-methylated *E.coli*. Therefore some restriction sites are blocked by methylation. If you wish to digest the vector using such sites you will need to transform the vector into a dam⁻ host and make fresh DNA.

Expression in mammalian cells

pFusionRed-N vector can be transfected into mammalian cells by any known transfection method. CMV promoter provides strong, constitutive expression of FusionRed or its fusions in eukaryotic cells. If required, stable transformants can be selected using G418 [Gorman 1985].

Propagation in E. coli

Suitable host strains for propagation in *E. coli* include DH5alpha, HB101, XL1-Blue, and other general purpose strains. Plasmid incompatibility group is pMB1/ColE1. The vector confers resistance to kanamycin (30 μ g/ml) to *E. coli* hosts. Copy number in *E. coli* is about 500.

Notice to Purchaser:

FusionRed-related materials (also referred to as "Products") are intended for research use only.

The Products are covered by European Pat. 1994149 and other Evrogen Patents and/or Patent applications pending. By use of these Products, you accept the terms and conditions of the applicable Limited Use Label License #001: http://www.evrogen.com/products/Evrogen-FP-license.shtml.

The CMV promoter is covered under U.S. Patents 5,168,062 and 5,385,839, and its use is permitted for research purposes only. Any other use of the CMV promoter requires a license from the University of Iowa Research Foundation, 214 Technology Innovation Center, Iowa City, IA 52242.